

www.mesoscale.com[®]

Ordering Information

MSD Customer Service
 Phone: 1-240-314-2795
 Fax: 1-301-990-2776
 Email: CustomerService@mesoscale.com

Scientific Support

Phone: 1-240-314-2798
 Email: ScientificSupport@mesoscale.com


Company Address

Meso Scale Discovery
 A division of
 Meso Scale Diagnostics, LLC.
 1601 Research Boulevard
 Rockville, MD 20850-3173 USA

Product Options	Catalog Number	Description
Multiplex	K15067M, K25067M	U-PLEX Biomarker Group 1 (human)
	K151AEM, K251AEM	U-PLEX Immuno-Oncology Group 1 (human)
	K151ACM, K251ACM	U-PLEX Metabolic Group 1 (human)
Singleplex	K151VHK-1/-2/-4	U-PLEX Human IFN- α 2a Assay with SECTOR™ plates
	K151VHK-21/-22/-24	U-PLEX Human IFN- α 2a Assay with QuickPlex Ultra™ plates
	K251VHK-2/-4	U-PLEX Human IFN- α 2a Assay with 384-well plates
Antibody Set	B21VH-2/-3	U-PLEX Human IFN- α 2a Antibody Set
Protocol	U-PLEX Product Inserts are available at http://www.mesoscale.com	

The MESO SCALE DISCOVERY[®] U-PLEX platform was designed to provide ultimate flexibility for detection of biomarkers in a wide variety of sample types. This datasheet provides the representative performance of the U-PLEX[®] Human IFN- α 2a Assay tested on U-PLEX 96-well SECTOR plates run as a multiplex. The data do not represent the product specifications. Under your experimental conditions, the assay may perform differently from the representative data. U-PLEX assays are offered in either singleplex or multiplex; both are available on 96- or 384-well plates. See a U PLEX product insert for instrument compatibility.

Representative Calibration Curve and Sensitivity

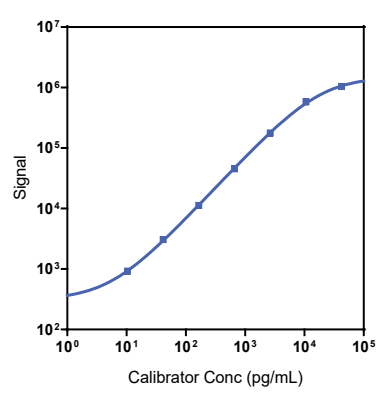

Assay	Median LLOD (pg/mL)	LLOD Range (pg/mL)
IFN- α 2a	4.0	1.9-6.4

The Calibrator curve was fitted with a 4-parameter logistic model with a $1/Y^2$ weighting. The lower limit of detection (LLOD) is a calculated concentration corresponding to 2.5 standard deviations above the background (zero Calibrator).

Precision

Control	Average Conc. (pg/mL)	Average Intra-run Conc. (%CV)	Inter-run Conc. (%CV)
High	11,500	9.9	9.0
Mid	1,090	7.1	9.7
Low	116	9.8	107

Controls were made by spiking Calibrator into assay diluent at 3 levels within the quantitative range of the assay. Average intra-run concentration %CV is the average %CV of the control replicates within an individual run. Inter-run concentration %CV is the variability of controls across multiple runs.

For Research Use Only.
 Not for use in diagnostic procedures.

MSD® U-PLEX Human IFN- α 2a

Tested Samples

Sample Type	Serum (N=10)	Plasma (N=10)	Spiked Plasma (N=5)	Spiked Serum (N=5)
Median (pg/mL)	NA	NA	292	277
Range (pg/mL)	NA	NA	266-337	248-321
% Detected	0	0	100	100

Normal serum and plasma samples were tested without dilution prior to the assay. NA = not applicable due to 0% detected

Dilution Linearity

Serum			EDTA Plasma		
Fold Dilution	Average % Recovery	% Recovery Range	Fold Dilution	Average % Recovery	% Recovery Range
2	120	106-152	2	113	103-119
4	119	104-163	4	113	100-119
8	133	109-181	8	111	88-120

Normal human serum and EDTA plasma were spiked with Calibrator and tested at different dilutions. Undiluted samples were tested to determine the expected concentration of the analyte. Samples may benefit from additional dilution with assay diluent to reduce matrix effects.

$$\% \text{ Recovery} = (\text{measured concentration} / \text{expected concentration}) \times 100$$

Spike Recovery

Spike Level	Serum		EDTA Plasma	
	Average % Recovery	% Recovery Range	Average % Recovery	% Recovery Range
High	82	74-98	96	83-108
Mid	83	76-92	84	74-89
Low	84	79-95	79	76-82

Normal serum and plasma were spiked with Calibrator at 3 levels. Undiluted samples were tested to determine the expected concentration of the analyte. Samples may benefit from additional dilution with assay diluent to reduce matrix effects.

$$\% \text{ Recovery} = (\text{measured concentration} / \text{expected concentration}) \times 100$$

Specificity

To assess specificity, the IFN- α 2a Antibody Set was tested individually against a larger panel of analytes for nonspecific binding (APRIL/TNFSF13, BAFF-R/TNFRSF13C, BCMA/TNFRSF17, BDNF, C-Peptide, CD20, CD27, CD28, CD40L (soluble), CD276/B7-H3, CTACK, CTLA-4, Desghrelin, ENA-78, Eotaxin, Eotaxin-2, Eotaxin-3, EPO, E-Selectin, FGF (basic), FGF-23, FLT3L, Fractalkine, FSH, Galectin-9, G-CSF, GITRL/TNFSF18, GTR/TNFRSF18, Ghrelin (Ser3-octanoylated), gp130 (soluble), GIP (1-42), GIP (3-42), GLP-1 (7-36), GLP-1 (9-36), GM-CSF, Granzyme A, Granzyme B, GRO- α , HAVCR2/TIM-3, HVEM/TNFRSF14, ICOS, ICOS-L/B7-H2, I-309, IFN- α 2a, IFN- β , IFN- γ , IL-1 α , IL-1 β , IL-1RA, IL-2, IL-2R α , IL-3, IL-4, IL-5, IL-6, IL-7, IL-8, IL-9, IL-10, IL-12/IL-23p40, IL-12p70, IL-13, IL-15, IL-16, IL-17A, IL-17A/F, IL-17C, IL-17D, IL-17E/IL-25, IL-17F, IL-18, IL-21, IL-22, IL-23, IL-27, IL-29/IFN- λ 1, IL-31, IL-33, Insulin, IP-10, LAG3, Leptin, LH, LIGHT/TNFSF14, MCP-1, MCP-2, MCP-4, M-CSF, MDC, MIF, MIG, MIP-1 α , MIP-1 β , MIP-5, MMP-1, MMP-2, MMP-7, Nectin-4, OX40/TNFRSF4, PD1, PD-L1, PD-L2, Pentraxin 3, Perforin, PIGF, PP, Proinsulin, proMMP-9, P-Selectin, PYY (3-36), RAGE (soluble), RANKL/TNFSF11, RANTES, S100A12, SDF-1 α , Tie-2, TIGIT, TLR1, TNF- α , TNF- β , TNF-RI, TNF-RII, TPO, TRAIL, TSLP, VEGF-A, VEGF-D, VEGFR-1/Flt-1, and YKL-40). Nonspecific binding was less than 2.0%.

$$\% \text{ Nonspecificity} = (\text{nonspecific signal} / \text{specific signal}) \times 100$$

Diluent Compatibility

Diluents 57 and 3 are provided with this assay. MSD offers a range of assay and antibody diluents for separate purchase. Depending on your assay needs, other diluents may be tested.

Assay Components

Calibrator: IFN- α 2a is included in Calibrator 3. The IFN- α 2a Calibrator is a full-length recombinant protein expressed in *E. coli*.

Antibodies: The U-PLEX Human IFN- α 2a Assay uses a mouse monoclonal antibody for capture and a mouse monoclonal antibody for detection.

Assay generation: A

Note: This datasheet contains representative assay performance data. In custom multiplex formats, the assay may perform differently from the representative data shown.

