

www.mesoscale.com<sup>®</sup>

**Ordering Information**

MSD Customer Service  
 Phone: 1-240-314-2795  
 Fax: 1-301-990-2776  
 Email: CustomerService@mesoscale.com

**Scientific Support**

Phone: 1-240-314-2798  
 Email: ScientificSupport@mesoscale.com


**Company Address**

Meso Scale Discovery  
 A division of  
 Meso Scale Diagnostics, LLC.  
 1601 Research Boulevard  
 Rockville, MD 20850-3173 USA

Product Options	Catalog Number	Description
<b>Multiplex</b>	K151AEM, K251AEM	U-PLEX Immuno-Oncology Group 1 (human)
<b>Singleplex</b>	K151AQQK-1/-2/-4	U-PLEX Human BTLA Assay with SECTOR™ plates
	K151AQQK-21/-22/-24	U-PLEX Human BTLA Assay with QuickPlex Ultra™ plates
	K251AQQK-2/-4	U-PLEX Human BTLA with 384-well plates
<b>Antibody Set</b>	B21AQQ-2/-3	U-PLEX Human BTLA Antibody Set
<b>Protocol</b>	U-PLEX Product Inserts are available at <a href="http://www.mesoscale.com">www.mesoscale.com</a>	

The Meso Scale Discovery U-PLEX platform was designed to provide ultimate flexibility for detection of biomarkers in a wide variety of sample types. This datasheet provides the representative performance of The MESO SCALE DISCOVERY<sup>®</sup> U-PLEX Human BTLA Assay tested on U-PLEX plates run as a multiplex. The data do not represent the product specifications. Under your experimental conditions, the assay may perform differently from the representative data. U-PLEX assays are offered in either singleplex or multiplex; both are available on 96- or 384-well plates. See a U-PLEX product insert for instrument compatibility.

**Representative Calibration Curve and Sensitivity**

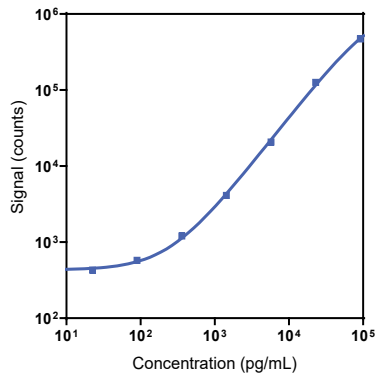

Assay	Median LLOD (pg/mL)	LLOD Range (pg/mL)
BTLA	60	48–132

The Calibrator curve was fitted with a 4-parameter logistic model with a 1/Y<sup>2</sup> weighting. The lower limit of detection (LLOD) is a calculated concentration corresponding to 2.5 standard deviations above the background (zero Calibrator).

**Precision**

Control	Average Conc. (pg/mL)	Average Intra-run Conc. (%CV)	Inter-run Conc. (%CV)
High	46,404	2.0	3.8
Mid	18,552	3.1	4.3
Low	7,533	3.2	5.0

Controls were made by spiking Calibrator into assay diluent at 3 levels within the quantitative range of the assay. Average intra-run concentration %CV is the average %CV of the control replicates within an individual run. Inter-run concentration %CV is the variability of controls across multiple runs.

For Research Use Only.  
 Not for use in diagnostic procedures.

# MSD® U-PLEX Human BTLA

## Tested Samples

Sample Type	Serum (N = 8)	EDTA Plasma (N = 8)	Citrate Plasma (N = 8)	Normal Lysate (N = 5)	Tumor Lysate (N = 15)
Median (pg/mL)	ND	513	451	256	482
Range (pg/mL)	42	ND–513	ND–451	ND–256	ND–995
% Detected	13	13	13	20	27

Normal serum and plasma samples were diluted 4-fold prior to the assay. Lysates were tested at a protein concentration of 0.5 mg/mL. ND = non-detectable (<LLOD).

## Dilution Linearity

Serum			EDTA Plasma		
Fold Dilution	Average % Recovery	% Recovery Range	Fold Dilution	Average % Recovery	% Recovery Range
2	106	76–118	2	113	85–124
8	87	64–96	8	90	74–99
16	70	52–83	16	76	65–88

Normal human serum and EDTA plasma were spiked with Calibrator and tested at different dilutions. Percent recovery at each dilution level was normalized to the dilution-adjusted, 4-fold concentration. Samples may benefit from additional dilution with assay diluent to reduce matrix effects.

$$\% \text{ Recovery} = (\text{measured concentration} / \text{expected concentration}) \times 100$$

## Spike Recovery

Spike Level	Serum		EDTA Plasma	
	Average % Recovery	% Recovery Range	Average % Recovery	% Recovery Range
High	90	76–95	111	95–119
Mid	93	80–98	105	91–111
Low	93	80–101	103	93–113

Normal serum and plasma were spiked with Calibrator at 3 levels. Spiked samples were diluted 4-fold to determine the expected concentration of the analyte. Samples may benefit from additional dilution with assay diluent to reduce matrix effects.

$$\% \text{ Recovery} = (\text{measured concentration} / \text{expected concentration}) \times 100$$

## Specificity

The BTLA Antibody Set was tested for nonspecific binding against all of the analytes in the Immuno-Oncology Group 1 and the majority of analytes in Biomarker Group 1. Any cross-reactivity greater than 2.0% is noted below. The U-PLEX Assay Designer shows compatible assays.

$$\% \text{ Nonspecificity} = (\text{nonspecific signal} / \text{specific signal}) \times 100$$

## Diluent Compatibility

Diluents 58 and 3 are provided when this is ordered in singleplex and multiplex assays.

## Assay Components

**Calibrator:** BTLA is included in Calibrator 32. The human BTLA Calibrator is BTLA (931–134) recombinant protein expressed in a human cell line.

**Antibodies:** The MESO SCALE DISCOVERY® U-PLEX Human BTLA Assay uses a mouse monoclonal antibody for capture and a mouse monoclonal antibody for detection.

**Assay generation:** A

**Note:** This datasheet contains representative assay performance data. In custom multiplex formats, the assay may perform differently from the representative data shown.

